

## NANO-MECHANICAL CHARACTERIZATION OF FREELY-SUSPENDED TWO-DIMENSIONAL MATERIALS WITH ATOMIC FORCE MICROSCOPY

M. Dimitropoulos<sup>1,2,\*</sup>, M. G. Pastore Carbone<sup>1</sup>, A. C. Manikas<sup>1,2</sup>, C. Galiotis<sup>1,2</sup>

<sup>1</sup> Institute of Chemical Engineering Sciences, Foundation for Research and Technology Hellas (FORTH/ICE-HT), 26504 Patras, Greece

<sup>2</sup> University of Patras, Chemical Engineering Department, 26504 Patras, Greece


(\*[mdim@chemeng.upatras.gr](mailto:mdim@chemeng.upatras.gr))

### ABSTRACT

The first isolation of graphene in 2004 by Novoselov and Geim<sup>[1]</sup>, has offered a new and exciting platform with extraordinary potential. The surge of two-dimensional (2D) materials has taken the world by storm and became a new research hotspot as the characterization of nanoscale properties became a necessity<sup>[2]</sup>. Lately, emphasis has been given to the synthesis of these materials and their heterostructures for mass production and integration in various applications, but the values of their properties are mostly theoretical and experimental data at the nanoscale are scarce. Therefore, in order to incorporate two-dimensional materials in the application of interest, one must first understand the full extent of their properties at the nanometer scale.

One of the most powerful techniques for nanoscale characterization is Atomic Force Microscopy (AFM), which can be practiced even at atomic scale resolution. The plethora of its sub-modes can identify mechanical, electrical, magnetic, chemical and even vibrational properties at the nanometer scale, in conjunction with topography. More specifically, for mechanical characterization, Quantitative nano-Mechanical characterization along with Nanoindentation allow for modulus/adhesion mapping and measurements, respectively. The most widely used approach for mechanical properties examination is obtained from force curves as a function of the indentation depth at the center of a freely-suspended membrane (**Figure**)<sup>[3]</sup>. The analysis models take into consideration the existing pre-tension and the bending rigidity of the membrane, as well as the Young's modulus and the hardness of the material<sup>[4]</sup>.

The purpose of this study is to investigate the mechanical properties of suspended 2D materials and their heterostructures on top of patterned substrates, with the aid of AFM. Initially, various 2D materials were directly mechanically exfoliated or transferred from CVD grown samples onto patterned SiO<sub>2</sub>/Si substrates. Subsequently, the flakes were identified for their quality and quantity with optical microscopy and Raman spectroscopy. Finally, the freely-suspended membranes of 2D materials allowed us to attain Young's modulus measurements without the contribution of the substrate, which were in close accordance with the theoretical predictions, as well as interface and layer-by-layer interactions of heterostructures. This method permits the evaluation of the theoretical data and can corroborate their significance.

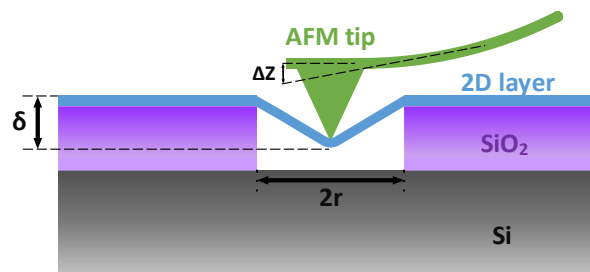

**Figure:** Typical nanoindentation experiment.

### REFERENCES

- [1] K.S. Novoselov, A.K. Geim, S. V. Morozov, Y.Z. D. Jiang, S. V. Dubonos, I. V. Grigorieva, A.A. Firsov, *Science* (80-. ). 306 (2004) 666–669.
- [2] K.S. Novoselov, A. Mishchenko, A. Carvalho, A.H. Castro Neto, *Science* (80-. ). 353 (2016).
- [3] A. Castellanos-Gomez, V. Singh, H.S.J. Van Der Zant, G.A. Steele, *Ann. Phys.* 527 (2015) 27–44.
- [4] C. Androulidakis, K. Zhang, M. Robertson, S.H. Tawfick, *2D Mater.* (2018).